DEC 2 9 2006

TITLE: CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF CARDIAC PERFORMANCE USING HEART RATE

INVENTOR'S NAME: RODNEY W. SALO, ET AL.

SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1

/0

REPLACEMENT SHEET

1/8

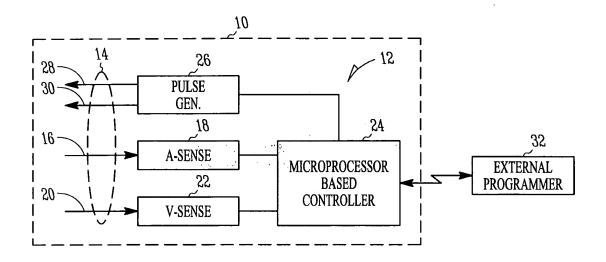


FIG. 1

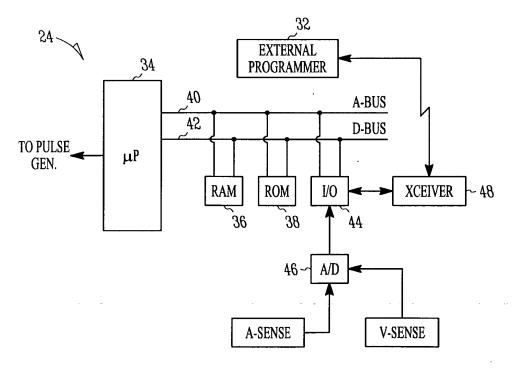


FIG. 2

TITLE: CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF CARDIAC PERFORMANCE USING HEART RATE

INVENTOR'S NAME: RODNEY W. SALO, ET AL.

SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1

2/8

REPLACEMENT SHEET

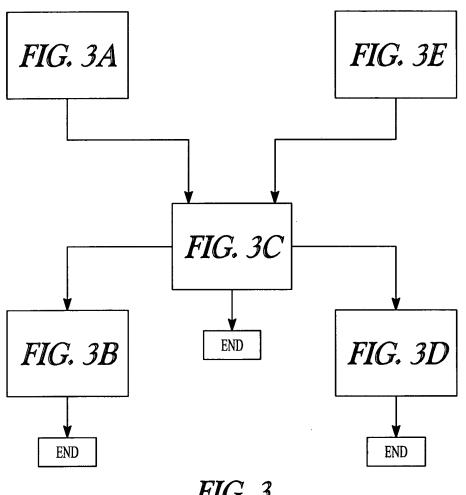


FIG. 3

TITLE: CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF

CARDIAC PERFORMANCE USING HEART RATE INVENTOR'S NAME: RODNEY W. SALO, ET AL.

SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1 REPLACEMENT SHEET

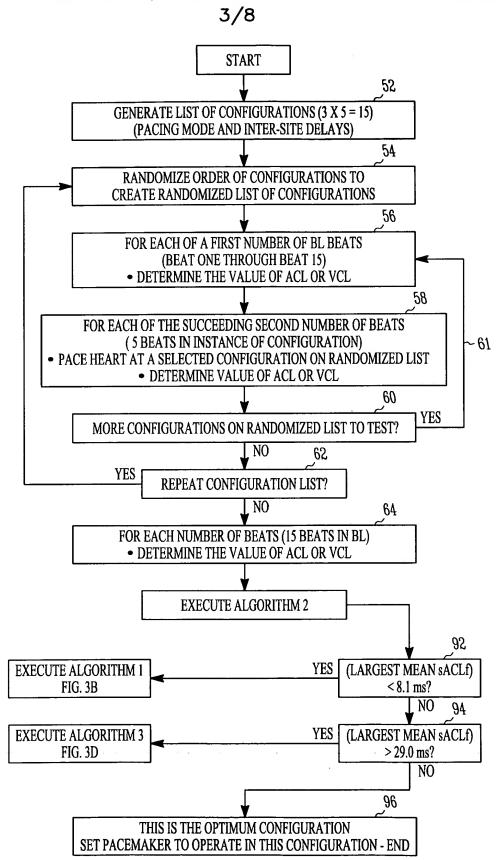


FIG. 3A

TITLE: CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF CARDIAC PERFORMANCE USING HEART RATE

INVENTOR'S NAME: RODNEY W. SALO, ET AL. SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1

4/8

REPLACEMENT SHEET

74

76

ALGORITHM 1

SMOOTH THE ARRAY OF ACLS OR VCLS WITH 3 POINT MOVING RECTANGLE WINDOW TO CREATE SMOOTHED ARRAY: SACL OR SVCL

FOR EACH CONFIGURATION INSTANCE:

IF ANY ABNORMAL BEATS (PVC, ETC.) FROM 8 BEATS BEFORE FIRST PACED BEAT TO 8 BEATS AFTER LAST PACED BEAT, THEN THIS IS DEFINED TO BE AN INVALID INSTANCE, ELSE THIS IS DEFINED TO BE A VALID INSTANCE

FOR EACH CONFIGURATION:

IF LESS THAN 3 VALID INSTANCES THEN THIS IS DEFINED TO BE AN INVALID CONFIGURATION, ELSE THIS IS DEFINED TO BE A VALID CONFIGURATION

68

FOR EACH VALID CONFIGURATION VALID INSTANCE:

- DETERMINE MAX. VALUE OF SACL OR SVCL IN INTERNAL FROM 2 BEATS AFTER FIRST BEAT OF CONFIGURATION INSTANCE TO 7 BEATS AFTER FIRST BEAT OF CONFIGURATION INSTANCE
- DETERMINE MEAN VALUE OF sACL OR SVCL IN INTERNAL FROM 2 BEATS BEFORE FIRST BEAT OF CONFIGURATION INSTANCE TO 2 BEATS AFTER FIRST BEAT OF CONFIGURATION INSTANCE
- DETERMINE SACL OR SVCL FEATURE VALUE SACLF OR SVCLF: IS MAGNITUDE (MAX.-MEAN) > MAGNITUDE (MIN.-MEAN) ? YES: SACLF OR SVCLF = (MAX.-MEAN) ~ 70 NO: SACLF OR SVCLF = (MIN.-MEAN) ~ 72

FOR EACH OF THE VALID CONFIGURATIONS:

• DETERMINE MEAN OF SACLE OR SYCLE OVER THE VALID INSTANCES
OF THE CONFIGURATION

AMONG VALID CONFIGURATIONS:

DETERMINE THE CONFIGURATION WITH LARGEST MEAN OF SACLE OR SVCLE
 THIS CONFIGURATION IS OPTIMUM CONFIGURATION.
 SET PACEMAKER TO OPERATE IN THIS OPTIMUM CONFIGURATION

FIG. 3B

SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1 REPLACEMENT SHEET 5/8 **ALGORITHM 2** 77 SMOOTH THE ARRAY OF ACLS OR VCLs WITH 11 POINT MOVING BLACKMAN WINDOW TO CREATE SMOOTHED ARRAY: SACL OR SVCL 78 FOR EACH CONFIGURATION INSTANCE: IF ANY ABNORMAL BEATS (PVC, ETC.) FROM 8 BEATS BEFORE FIRST PACED BEAT TO 8 BEATS AFTER LAST PACED BEAT, THEN THIS IS DEFINED TO BE AN INVALID INSTANCE, ELSE THIS IS DEFINED TO BE A VALID INSTANCE 79 FOR EACH CONFIGURATION: IF LESS THAN 3 VALID INSTANCES THEN THIS IS DEFINED TO BE AN INVALID CONFIGURATION, ELSE THIS IS DEFINED TO BE A VALID CONFIGURATION 90 FOR EACH VALID CONFIGURATION VALID INSTANCE: DETERMINE MAX. VALUE OF SACL OR SVCL IN INTERNAL FROM 2 BEATS AFTER FIRST BEAT OF CONFIGURATION INSTANCE TO 8 BEATS AFTER FIRST BEAT OF CONFIGURATION INSTANCE DETERMINE MIN. VALUE OF SACL OR SVCL IN INTERNAL FROM 3 BEATS BEFORE FIRST BEAT OF CONFIGURATION INSTANCE TO BEAT WITH MAX. VALUE DETERMINE SACL OR SVCL FEATURE VALUE SACLE OR SVCLE: SACLE OR SVCLE = MAX. -MIN. 85 FOR EACH OF THE VALID CONFIGURATIONS: DETERMINE MEAN OF SACLE OR SVCLE OVER THE VALID INSTANCES 84 AMONG THE VALID CONFIGURATIONS: DETERMINE CONFIGURATION WITH THE LARGEST MEAN SACLF OR SVCLF 98 AMONG THE VALID INSTANCES OF THIS VALID CONFIGURATION: DETERMINE MEDIAN VALUE OF sACLf OR sVCLf YES DETERMINE MAX. VALUE OF sACLF OR sVCLF 91 -IS (MAX./MEDIAN-1) < PREDEFINED THRESHOLD? 90 88 THE INSTANCE WITH THE MAX. VALUE OF SACLE OR SVCLE IS DEFINED TO BE THIS IS THE AN INVALID INSTANCE. IF THIS VALID CONFIGURATION NOW HAS < 3 VALID CANDIDATE OPTIMUM INSTANCES, THEN IT IS DEFINED TO BE AN INVALID CONFIGURATION. CONFIGURATION ELSE CALCULATE MEAN OF VALUES OF VALID INSTANCES OF THIS VALID CONFIGURATION

TITLE: CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF CARDIAC PERFORMANCE USING HEART RATE
. INVENTOR'S NAME: RODNEY W. SALO, ET AL.

FIG. 3C

CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF CARDIAC PERFORMANCE USING HEART RATE INVENTOR'S NAME: RODNEY W. SALO, ET AL. SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1 REPLACEMENT SHEET 6/8 **ALGORITHM 3** SMOOTH THE ARRAY OF ACLS OR VCLS WITH 11 POINT MOVING BLACKMAN WINDOW TO CREATE SMOOTHED ARRAY: SACL OR SVCL FOR EACH CONFIGURATION INSTANCE: IF ANY ABNORMAL BEATS (PVC, ETC.) FROM 8 BEATS BEFORE FIRST PACED BEAT TO 8 BEATS AFTER LAST PACED BEAT, THEN THIS IS DEFINED TO BE AN INVALID INSTANCE, ELSE THIS IS DEFINED TO BE A VALID INSTANCE FOR EACH CONFIGURATION: IF LESS THAN 3 VALID INSTANCES THEN THIS IS DEFINED TO BE AN INVALID CONFIGURATION, ELSE THIS IS DEFINED TO BE A VALID CONFIGURATION 97 FOR EACH VALID CONFIGURATION VALID INSTANCE: DETERMINE MAX. VALUE OF SACL OR SVCL IN INTERNAL FROM 1 BEAT AFTER FIRST BEAT OF CONFIGURATION INSTANCE TO 8 BEATS AFTER FIRST BEAT OF CONFIGURATION INSTANCE DETERMINE MIN. VALUE OF SACL OR SVCL IN INTERNAL FROM 3 BEATS BEFORE FIRST BEAT OF CONFIGURATION INSTANCE TO BEAT WITH MAX. VALUE • DETERMINE SACL OR SVCL FEATURE VALUE SACLE OR SVCLE. sACLf OR sVCLf = MAX.-MIN. 98 FOR EACH OF THE VALID CONFIGURATIONS: DETERMINE MEAN OF SACLE OR SVCLE OVER THE VALID INSTANCES AMONG THE VALID INSTANCES OF THIS VALID CONFIGURATION: DETERMINE MEDIAN VALUE OF sACLf OR sVCLf NO. DETERMINE MAX. VALUE OF sACLI OR sVCLI


IS (MAX./MEDIAN-1) < PREDEFINED THRESHOLD? YES 100 THIS IS THE OPTIMUM CONFIGURATION. SET PACEMAKER TO OPERATE IN THIS CONFIGURATION

FIG. 3D

TITLE: CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF

CARDIAC PERFORMANCE USING HEART RATE INVENTOR'S NAME: RODNEY W. SALO, ET AL.

SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1 REPLACEMENT SHEET

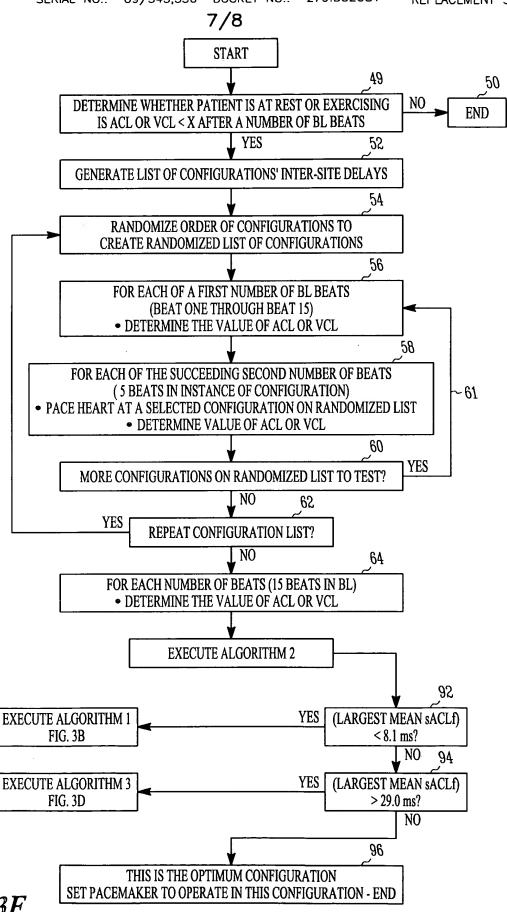

CARDIAC RHYTHM MANAGEMENT SYSTEM WITH OPTIMIZATION OF CARDIAC PERFORMANCE USING HEART RATE

INVENTOR'S NAME: RODNEY W. SALO, ET AL.

SERIAL NO.: 09/545,536 DOCKET NO.: 279.D32US1 REPLACEMENT SHEET

FIG. 5

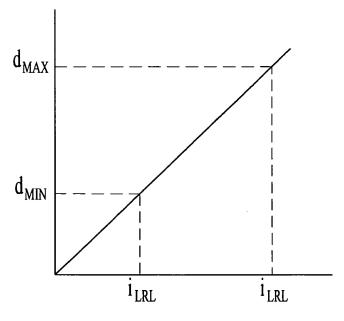


FIG. 6